

# Resolving the Complexity of Spatial Lipidomics Using MALDI TIMS Imaging Mass Spectrometry

Katerina Djambazova, Dustin R. Klein, [Lukasz Migas](#), Elizabeth Neumann, Emilio Rivera, Raf Van de Plas, Richard M. Caprioli, [Jeffrey Spraggins](#)

Submitted date: 13/06/2020 · Posted date: 16/06/2020

Licence: CC BY-NC-ND 4.0

Citation information: Djambazova, Katerina; Klein, Dustin R.; Migas, Lukasz; Neumann, Elizabeth; Rivera, Emilio; Van de Plas, Raf; et al. (2020): Resolving the Complexity of Spatial Lipidomics Using MALDI TIMS Imaging Mass Spectrometry. ChemRxiv. Preprint. <https://doi.org/10.26434/chemrxiv.12331652.v2>

Lipids are a structurally diverse class of molecules with important biological functions including cellular signaling and energy storage. Matrix-assisted laser desorption/ionization (MALDI) imaging mass spectrometry (IMS) allows for direct map-ping of biomolecules in tissue. Fully characterizing the structural diversity of lipids remains a challenge due to the presence of isobaric and isomeric species, which greatly complicates data interpretation when only  $m/z$  information is available. Integrating ion mobility separations aids in deconvoluting these complex mixtures and addressing the challenges of lipid IMS. Here we demonstrate that a MALDI quadrupole time-of-flight (Q-TOF) mass spectrometer with trapped ion mobility spectrometry (TIMS) enables approximately a ~270% increase in the peak capacity during IMS experiments. MALDI TIMS-MS separation of lipid isomer standards, including sn-backbone isomers, acyl chain isomers, as well as double bond positional and geometric isomers are demonstrated. As a proof-of-concept, in situ separation and imaging of lipid isomers with distinct spatial distributions was performed using tissue sections from a whole-body mouse pup.

## File list (2)

---

0612\_Manuscript\_v2.pdf (0.96 MiB)

[view on ChemRxiv](#) • [download file](#)

---

0612\_Supplemental.pdf (889.59 KiB)

[view on ChemRxiv](#) • [download file](#)

---

# Resolving the Complexity of Spatial Lipidomics Using MALDI TIMS Imaging Mass Spectrometry

Katerina V. Djambazova<sup>1,2</sup>, Dustin R. Klein<sup>2,3</sup>, Lukasz G. Migas<sup>4</sup>, Elizabeth K. Neumann<sup>2,3</sup>, Emilio S. Rivera<sup>2,3</sup>, Raf Van de Plas<sup>4</sup>, Richard M. Caprioli<sup>1-3,5,6</sup> and Jeffrey M. Spraggins<sup>1-3\*</sup>.

<sup>1</sup>Department of Chemistry, Vanderbilt University, 7330 Stevenson Center, Station B 351822, Nashville, TN 37235, USA

<sup>2</sup>Mass Spectrometry Research Center, Vanderbilt University, 465 21<sup>st</sup> Ave S #9160, Nashville, TN 37235, USA

<sup>3</sup>Department of Biochemistry, Vanderbilt University, 607 Light Hall, Nashville, TN 37205, USA

<sup>4</sup>Delft Center for Systems and Control, Delft University of Technology, 2628 CD Delft, The Netherlands

<sup>5</sup>Department of Pharmacology, Vanderbilt University, 2220 Pierce Avenue, Nashville, TN 37232, USA

<sup>6</sup>Department of Medicine, Vanderbilt University, 465 21<sup>st</sup> Ave S #9160, Nashville, TN 37235, USA

**KEYWORDS.** *Imaging mass spectrometry, Matrix-assisted laser desorption/ionization, Trapped ion mobility spectrometry, Ion mobility mass spectrometry, Lipidomics*

---

**ABSTRACT:** Lipids are a structurally diverse class of molecules with important biological functions including cellular signaling and energy storage. Matrix-assisted laser desorption/ionization (MALDI) imaging mass spectrometry (IMS) allows for direct mapping of biomolecules in tissue. Fully characterizing the structural diversity of lipids remains a challenge due to the presence of isobaric and isomeric species, which greatly complicates data interpretation when only  $m/z$  information is available. Integrating ion mobility separations aids in deconvoluting these complex mixtures and addressing the challenges of lipid IMS. Here we demonstrate that a MALDI quadrupole time-of-flight (Q-TOF) mass spectrometer with trapped ion mobility spectrometry (TIMS) enables approximately a ~270% increase in the peak capacity during IMS experiments. MALDI TIMS-MS separation of lipid isomer standards, including sn-backbone isomers, acyl chain isomers, as well as double bond positional and geometric isomers are demonstrated. As a proof-of-concept, *in situ* separation and imaging of lipid isomers with distinct spatial distributions was performed using tissue sections from a whole-body mouse pup.

---

Lipids are crucial components of cellular membranes with essential functions in cellular structure, signaling, energy storage, and homeostasis.<sup>1-4</sup> They are structurally diverse molecules, with many having the same nominal mass and some the same exact mass. Developing more comprehensive ways of characterizing these complex mixtures in tissue is important since even subtle structural differences can have significant biological implications. Mass spectrometry (MS) is a powerful tool for lipid analysis with common *ex vivo* approaches utilizing electrospray ionization (ESI) and either direct infusion of lipid extracts or separations via liquid chromatography, neither of which retain spatial information of the analyte in tissue.<sup>3,5</sup> In contrast, imaging mass spectrometry (IMS), which analyzes biomolecules *in situ*, enables highly specific mass analysis while retaining analyte spatial information.<sup>6</sup> Matrix-assisted laser desorption/ionization (MALDI), one of the most commonly used ionization techniques for IMS, permits visualization of lipids<sup>7,8</sup>, metabolites<sup>9,10</sup>, glycans<sup>11,12</sup>, peptides<sup>14</sup>, and proteins<sup>15,16</sup>. In a typical MALDI IMS experiment, tissues are thinly sectioned, mounted on a conductive glass slide, and coated with an energy absorbing matrix. Desorption and ionization is facilitated by laser ablation at specific positions on the tissue surface. An individual mass spectrum is generated at each position (pixel), and molecular images are generated by plotting the ion intensity as

a heat map across all the pixels sampled on the tissue.<sup>16</sup> The spatial resolution of these chemical images is defined by the laser ablation area on the tissue and the spacing between the pixels (pitch).

The structural diversity of lipids in tissues can often be underrepresented during MALDI IMS experiments owing to the presence of isobars and isomers that cannot be distinguished by their measured  $m/z$  alone. High resolving power and high mass accuracy instrumentation, including Fourier transform ion cyclotron resonance (FT-ICR) and Orbitrap mass analyzers, have been critical in enabling both lipid identification at the sum composition level, as defined by Liebis et al<sup>2</sup>, and isobar differentiation.<sup>17,18</sup> However, since structural isomers have identical masses, they cannot be separated simply by high mass resolving power instruments. In these cases, ion images of single  $m/z$  values may depict composites of multiple lipid species. While tandem MS techniques (MS<sup>n</sup>) have been used to differentiate some isomers *in situ*, most methods using MS/MS are highly targeted and time-consuming.<sup>19</sup> Additionally, while chromatographic techniques are capable of separating isomeric lipid species<sup>20</sup>, they typically require lengthy acquisition times and are incompatible with the time-scales of IMS.

Ion mobility enables gas-phase separation of molecules according to their shape, size and charge on timescales ( $\mu\text{s}$ - $\text{ms}$ ) that are compatible with typical IMS experiments. Ion mobility separation techniques include drift tube ion mobility spectrometry (DTIMS)<sup>21,22</sup>, traveling wave ion mobility spectrometry (TWIMS)<sup>23,24</sup>, structures for lossless ion manipulations (SLIM)<sup>25,26</sup>, field asymmetric ion mobility spectrometry (FAIMS)<sup>27-29</sup>, and trapped ion mobility spectrometry (TIMS)<sup>30,31</sup>. DTIMS and TWIMS are operated under lower pressure regimes and have typically been coupled to time-of-flight (TOF) mass spectrometers. Mobility trendlines generated during DTIMS and TWIMS experiments have been beneficial for identifying classes and subclasses of lipids, and have also enabled increases in signal-to-noise measurements.<sup>22,23</sup> DTIMS and TWIMS have both been used to separate ESI-generated lipid isomers, including *sn*-positional, double bond positional, and stereochemical isomers (R versus S).<sup>32-34</sup> Furthermore, both TWIMS and DTIMS have previously been utilized on instruments capable of imaging mass spectrometry.<sup>35,36</sup> Recent developments in TWIMS technologies has led to the development of SLIM devices that are capable of high-resolution separations of many lipid isomers.<sup>38</sup> Although SLIM shows promise for elucidating subtle structural differences, it has only been coupled to liquid surface probes and has not yet been integrated with a MALDI IMS setup.<sup>40</sup> FAIMS devices do not require low pressures and can be integrated in a variety of mass spectrometers.<sup>41</sup> With the exception of a higher resolution planar-FAIMS device<sup>39</sup>, most FAIMS devices operate at relatively low resolving power and are used to filter specific classes or subclasses of lipids.<sup>33</sup> For example, shotgun lipidomic studies with FAIMS have been successful in lipid class separations<sup>33,42</sup>, and Lintonen et al. have shown the use of FAIMS for the construction of lipid class trendlines.<sup>43</sup> Similar to DTIMS and TWIMS, FAIMS has also been coupled to IMS platforms. For example, Feider et al. have demonstrated improved detection and imaging of cardiolipins in rat brain with the use of FAIMS and desorption electrospray ionization.<sup>44</sup>

Recently, TIMS integration with IMS has shown great potential for high-resolution separations.<sup>45</sup> TIMS platforms are capable of resolving powers  $>200$  in  $50 - 500$  ms.<sup>31,46-48</sup> Briefly, TIMS is performed in the first vacuum stage of the mass spectrometer in an augmented ion funnel. An electric field gradient (EFG) is applied to the ion tunnel to trap ions against a carrier gas and separate them based on mobility. Following ion trapping, the EFG is gradually reduced to sequentially elute ions with ascending mobilities.<sup>30,47</sup> The voltage range ( $\Delta V_{\text{ramp}}$ ) and ramp time ( $t_{\text{ramp}}$ ) dictate the scan rate ( $S_r = \Delta V_{\text{ramp}}/t_{\text{ramp}}$ )<sup>49</sup>, with slower scan rates resulting in higher resolution separations. While a MALDI TIMS-MS platform has successfully been used to separate and spatially map lipid species that overlap in  $m/z$ -space, a full assessment of the effect of TIMS on peak capacity and isomer separation in an imaging context has not yet been demonstrated.<sup>45,50</sup> In addition, the ability of TIMS to resolve isomeric lipids via liquid chromatography-TIMS-MS suggest it has the capability to resolve isomeric lipids *in situ*.<sup>50</sup>

Here we present the use of the TIMS technology for increased peak capacity and deconvolution of complex lipid mixtures in an imaging context. The high-resolution separation of MALDI-generated lipid isomer standards is demonstrated, followed by an *in situ* assessment of lipid peak capacity using TIMS-IMS and visualization of lipid isomers in full-body mouse pup tissue.

These data provide the foundation for maximizing sensitivity and specificity of lipid imaging experiments using TIMS-based MALDI platforms.

## METHODS

**Materials:** 1,5-diaminonaphthalene (DAN) and 2',5'-dihydroxyacetophenone (DHA) were purchased from Sigma-Aldrich (St. Louis, MO, USA). HPLC-grade acetonitrile, methanol, ethanol, chloroform, and tetrahydrofuran (THF) were purchased from Fisher Scientific (Pittsburgh, PA, USA). Lipid standards were purchased from Avanti Polar Lipids (Alabaster, AL, US).

**Sample Preparation:** For tissue preparation, one-week old C57BL/6 control mouse pup was snap frozen at  $-80$  °C, shaved over dry ice, and cryosectioned ( $20\mu\text{m}$  thickness), using a CM3050 S cryostat (Leica Biosystems, Wetzlar, Germany). The tissue was thaw-mounted onto conductive indium tin oxide coated glass slides (Delta Technologies, Loveland, CO, USA). The sample was placed on a heated metal block ( $75$  °C) and coated with a  $15$  mg/mL solution of DAN dissolved in THF using a TM Sprayer (HTX Technologies, LLC, Chapel Hill, NC, USA) ( $0.05$  mL/hr, 5 passes, 0 sec drying time,  $40$  °C spray nozzle). For lipid standard analysis, the standards (in chloroform) were aliquoted in vials ( $5\text{mg/mL}$ ), dried down with nitrogen gas, and dissolved in  $80\%$  methanol for a final concentration of  $1.0$  mg/mL. The aliquots were then mixed with DHA matrix, dissolved in  $90\%$  acetonitrile to a final lipid concentration of  $0.2$  mg/mL. The mixtures were spotted on a polished steel target – MTP AnchorChip (Bruker Daltonik, Bremen, Germany).

**MALDI TIMS-MS:** All experiments were carried out on a prototype timsTOF fleX mass spectrometer (Bruker Daltonik, Bremen, Germany).<sup>45</sup> Images were acquired in TIMS mode of operation with an ion transfer time of  $100$   $\mu\text{s}$ , pre-pulse storage time of  $8$   $\mu\text{s}$ , a collision RF of  $2000$  Vpp, a TIMS funnel 1 (accumulation) RF of  $450$  Vpp, a TIMS funnel 2 RF (analysis) of  $400$  Vpp, a multipole RF of  $400$  Vpp, and a collision cell entrance (in) voltage of  $300$  V. The data ( $196514$  pixels) were collected at  $50$   $\mu\text{m}$  spatial resolution with beam scan on ( $46$   $\mu\text{m}$ ), using  $200$  shots per pixel and  $48\%$  laser power. Data were collected in positive ion mode from  $m/z$   $50 - 1500$ . The TIMS EFG scan time was set to  $650$  ms, with reduced mobility range of  $1.370 - 1.520$  Vs/cm<sup>2</sup> ( $S_r = 0.03$  V/ms). TIMS imaging data was visualized using custom in-house developed software. Lipid identifications were determined with high mass accuracy and LIPIDMAPS lipidomics gateway (lipidmaps.org<sup>51,52</sup>).

**Standards:** Data were acquired in TIMS mode of operation with specific tuning of TIMS ramp time, mobility range, and scan rate for each of the analyzed standards (**Table S1**). Data were acquired at  $30$   $\mu\text{m}$  spatial resolution with  $\sim 75\%$  laser power at  $10$  kHz,  $200$  shots per pixel, and  $\sim 300$  pixels per sample. For all isomer standard experiments, the scan time was kept constant ( $650$  ms), and the mobility range was manipulated to control the voltage change as well as the scan rate.

**MALDI FT-ICR MS:** Whole-body mouse pup tissue data were acquired using a 15T MALDI FT-ICR mass spectrometer (Bruker Daltonics, Billerica, MA, USA). The images were generated with the “small” laser setting  $\sim 50$   $\mu\text{m}$ , with a spatial resolution of  $150$   $\mu\text{m}$ . The data was collected with  $300$  laser shots


per pixels at 46% laser power. The recorded mass spectrum had a  $m/z$  range of 460–1500 with a transient length of 4.2 s, resulting in  $\sim 600\,000$  resolving power at  $m/z$  760.

**Data Processing:** The mouse pup data was exported into a custom binary format optimized for storage and speed of analysis of the ion mobility-IMS data collection. Each frame/pixel contains between 10,000 – 100,000 centroid peaks that span the acquisition range of  $m/z$  50 – 1500 and  $1/K_0$  1.370–1.520  $\text{Vs}/\text{cm}^2$  with 400,577 and 5857 bins in the MS and ion mobility-dimensions, respectively. The processing pipeline requires common  $m/z$  and  $1/K_0$  axes, hence individual centroid peaks were inserted at their correct bin position along the MS and ion mobility-dimensions; missing values were set to zero. Following the conversion process, a mean mass spectrum of the entire dataset was generated and peak-picked. A total of  $\sim 400$  most intense ions were then selected and extracted to generate ion mobility rich ion images. Each of the extracted ion mobility-ion images was peak-picked (in the ion mobility-dimension) and automatically fitted with Gaussian-distributions to identify single and multi-conformational species. Ion mobility-selected ion images were then visualized to examine the conformational-specific localization in the spatial domain.

## RESULTS AND DISCUSSION

**TIMS Peak Capacity:** Direct sampling of lipids from tissue by imaging mass spectrometry produces complex mass spectra. To assess the advantage of incorporating TIMS into an IMS workflow, spectral peak capacity (i.e. the number of detectable features) was determined for an average spectrum of a whole-body mouse pup tissue image. Peak capacity is determined by peak picking the data to identify the number of unique mass spectral features with and without incorporating TIMS data. When the TIMS dimension is incorporated, the peak capacity is determined as the number of unique mass spectral and associated mobility features. The average mass spectrum from  $m/z$  680–880 with no ion mobility separation shown in **Figure 1A** has a peak capacity of 247. Over the same  $m/z$  range, with a 650 ms TIMS separation, the peak capacity increases to 916. **Figure 1B** provides a heat map representation of the ion mobility-MS data. Integrating ion mobility into the MALDI IMS workflow, therefore, provides a 271% increase in the number of species detected for lipid IMS experiments.


The increased peak capacity afforded by TIMS results from the separation of isobaric and isomeric species that otherwise overlap in the  $m/z$  dimension. A well-documented example of isobaric interference that occurs with lipid analysis results from overlap of an M+2 isotopologue of one species with the monoisotopic (M) peak of another lipid with one less double bond. This ‘double bond ambiguity’ gives rise to a mass difference of 0.00895 Da.<sup>53</sup> Peak overlap that results from unresolved isotopic distributions can lead to artificially high mass errors that preclude lipid identification.<sup>53,54</sup> While TOF instruments generally are not able to resolve peaks resulting from lipid double bond ambiguity, the added dimension of ion mobility allows for gas phase separations that help address this spectral complexity. **Figure 2** provides an example of how TIMS can resolve overlapping isotopes for the ion of  $m/z$  848.55. With a mobility separation time of 650 ms, the ion of  $m/z$  848.55 resolves into two mobility peaks at  $1/K_0$  1.485 – 1.499 and  $1/K_0$  1.468 – 1.485. **Figure 2A** shows the composite ion image combining the mobility separated species for  $m/z$  848.55. **Figures 2B** and **2C**

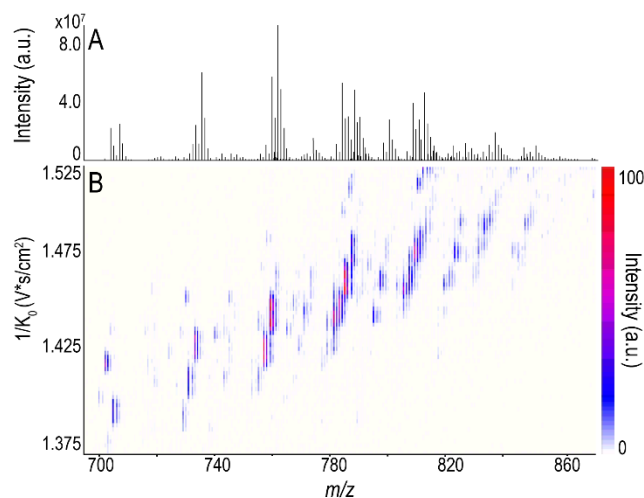

**Figure 1:** Average mass spectrum (**A**) and ion mobility heat map (**B**) generated for average whole-body mouse pup data highlight the presence of multiple isobaric and isomeric species. The diversity of ions observed in the heat map for each  $m/z$  highlight the complexity of the data and demonstrate the increased peak capacity by incorporating the TIMS separation dimension. The mass spectrum ranges include  $m/z$  values 680 – 880 and the heat map incorporates  $1/K_0$  values 1.375 – 1.500.

show the individual ion images for  $1/K_0$  1.485 – 1.499 and  $1/K_0$  1.468 – 1.485, respectively, each of which localizes to distinct regions of the tissue. For example,  $1/K_0$  1.485 – 1.499 is found at higher abundance in the brain and with more distinct localization to the intestinal walls as compared to  $1/K_0$  1.468 – 1.485, which shows unique localization in the adipose tissue. **Figure 2D** shows the ion image of  $m/z$  846.54 for the mobility range of  $1/K_0$  1.468 – 1.485, which, by comparison, localizes to the same regions of tissue as  $m/z$  848.55 for the mobility range of  $1/K_0$  1.468 – 1.485. Considering the intact mass measurements and differences in spatial localizations, the species at  $m/z$  848.55 with a mobility elution time of 1.485 – 1.499 likely corresponds to the monoisotopic peak of phosphatidylcholine (PC)  $[\text{PC}(38:4)+\text{K}]^+$ , while the species of  $m/z$  848.55 with a mobility elution time of 1.468 – 1.485 corresponds to the M+2 peak of  $[\text{PC}(38:5)+\text{K}]^+$ . Without a mobility separation,  $m/z$  848.55 would likely be assigned solely as the higher intensity  $[\text{PC}(38:4)+\text{K}]^+$  and the generated ion image would provide a skewed representation of its distribution throughout the tissue.

**Separation of MALDI-generated Lipid Isomer Standards:** Isomers have the same chemical formulas (identical  $m/z$ -values), and therefore cannot be separated by mass resolving power. Gas-phase separations prior to MS analysis have been successful in resolving isomeric species. Here MALDI TIMS isomer separations were evaluated in both positive and negative ion mode. The performance was assessed using standards including isomers with differences in sn-position, acyl chain composition, double bond position, and double bond geometry (**Figure S1**). These standards have previously been characterized using ESI techniques in combination with TIMS, as well as other separation techniques.<sup>34,55</sup> MALDI TIMS analyses, however, lay the groundwork for isomer differentiation within a molecular imaging context. The mass spectra of the individual standards can be found in **Figures S2** and **S3**. Extracted ion

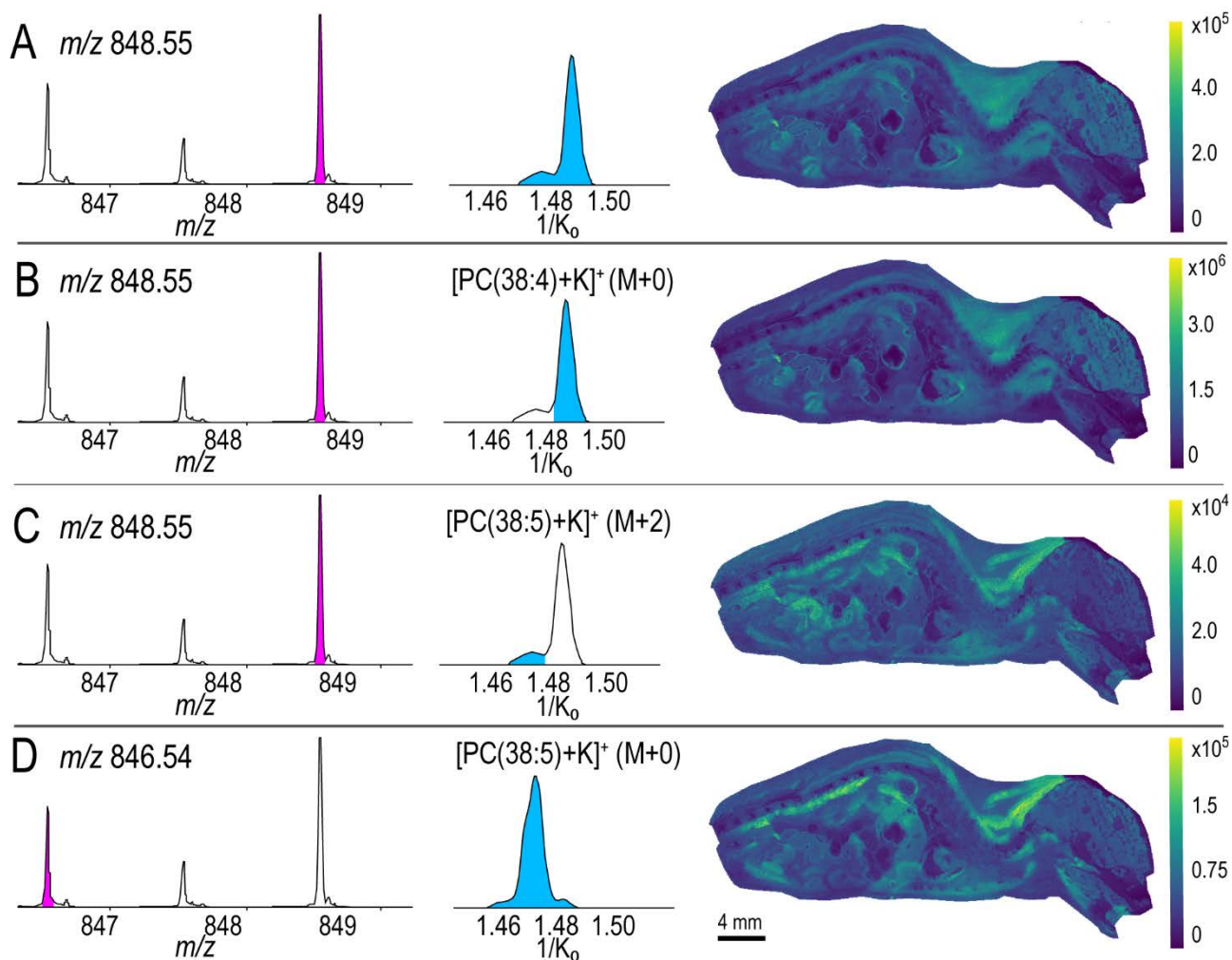

**Figure 2:** The added dimension of ion mobility helps address spectral ambiguities to show how the double bond ambiguity is resolved for the ion of  $m/z$  848.55. The ion of  $m/z$  848.55 resolves into two peaks at  $1/K_0$  1.485 – 1.499 and  $1/K_0$  1.468 – 1.485 (A). A composite image of the two peaks (A) and individual ion images for  $1/K_0$  1.485 – 1.499 (B) and  $1/K_0$  1.468 – 1.485 (C) highlight differences in their spatial distributions. The ion image of  $m/z$  846.54 ( $1/K_0$  1.468 – 1.485) (D) localizes to the same regions as  $m/z$  848.55 for the mobility range of  $1/K_0$  1.468 – 1.485 (C). The species of  $m/z$  848.55 with  $1/K_0$  1.485 – 1.499, likely corresponds to the M+0 (monoisotopic) peak of  $[PC(38:4)+K]^+$ , while the species of  $m/z$  848.55 with  $1/K_0$  1.468 – 1.485 corresponds to the M+2 peak of  $[PC(38:5)+K]^+$ .


mobilograms were generated for each standard as well as for mixtures of isomeric standards. Both ion mobility resolving power (R), and resolution (r), as defined by Fernandez-Lima and co-workers<sup>55</sup> were calculated where appropriate.

#### Sn-position

Subtle differences in glycerophospholipid structure, including differences in acyl chain configuration, can have dramatic effects on biological function. A unique case of sn-positional isomers is bis(monoacylglycerol)phosphates (BMPs) and phosphatidylglycerols (PGs). In contrast to PGs, where both acyl chains are connected to the glycerol backbone, BMPs have one acyl chain connected at the sn-1 position on the glycerol backbone and a second acyl chain connected to the sn-1' position of the glycerol headgroup, resulting in an unusual sn-1:sn-1' configuration.<sup>31,46</sup> Differentiation of BMPs and PGs is important, as both are present in mammalian cells, but differ in their biological roles and cell locations. BMP lipids have been found in the endosomal and lysosomal compartments of the cell and have

been associated with Niemann-Pick disease, an endosomal/lysosomal storage disease<sup>56,57</sup>. PGs, on the other hand, are thought to have antiviral properties, inhibiting viral proliferation in the lungs including respiratory syncytial virus.<sup>58,59</sup>

Sn-positional isomers are particularly challenging to characterize as collision-based tandem mass spectrometry methods produce spectra with only differences in ion abundances, i.e. non-unique identifying ions.<sup>59</sup> To evaluate MALDI TIMS as a platform to separate PG and BMP isomers, mobility data was collected for the standards PG (14:0/14:0) and BMP (14:0/14:0) (Figure 3A). Standards were first analyzed individually to determine mobility values. The MS1 spectrum for each standard contains a predominant ion of  $m/z$  665.49, corresponding to the  $[M-H]^-$  species for each lipid. The measured  $1/K_0$  values were 1.290 and 1.299 for PG (14:0/14:0) and BMP (14:0/14:0), respectively. The difference in backbone geometry results in an increased distance between the acyl chains and consequently an

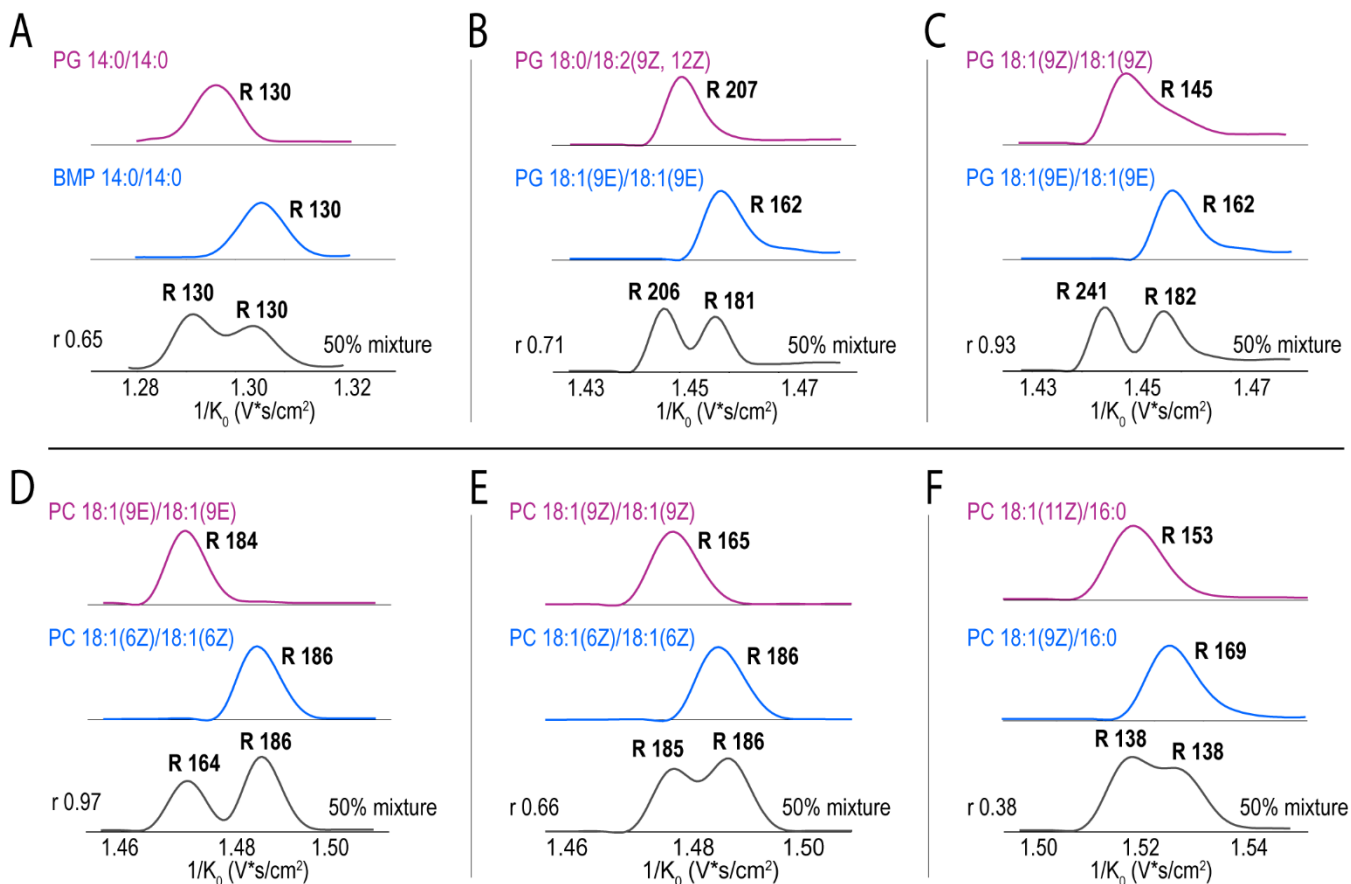

**Figure 3:** Extracted ion mobilograms of MALDI-generated ions of isomeric lipid standards are shown for both negative (A-C) and positive ionization mode (D-F). Negative ionization mode analysis of sn-backbone isomers PG (14:0/14:0) and BMP (14:0/14:0) ( $[M-H]^-$ ,  $m/z$  665.49) (A); Analysis of double bond position and geometry isomers – PG (18:0/18:2 (9Z, 12Z)) and PG (18:1/18:1 (9E)) (B) and PG (18:1(9Z)/18:1 (9Z)) and PG ((18:1(9E)/18:1(9E)) ( $[M-H]^-$ ,  $m/z$  773.53) (C). Positive ion mode analysis of double bond position and geometry isomers – PC (18:1(9E)/18:1(9E)) and PC (18:1(6Z)/18:1(6Z)) (D) and PC (18:1(9Z)/18:1(9Z)) and PC(18:1(6Z)/18:1(6Z)) ( $[M+Na]^+$ ,  $m/z$  808.57) (E), as well as PC (18:1(11Z)/16:0) and PC18:1(9Z)/16:0) ( $[M+Na]^+$ ,  $m/z$  782.56) (F). Resolving power (R) is listed where appropriate.

increased difference in collisional cross section between BMP and PG.<sup>34</sup> This is consistent with previous findings showing that PGs have more compact structures than the BMPs.<sup>34,59</sup> A 1:1 mixture was analyzed to evaluate the separation of the isomers. An extracted ion mobilogram was generated for the deprotonated species ( $m/z$  665.49), and two peaks could be resolved in the mobilogram (1.290 and 1.299). The mobilogram of the mixture shows separation of the two species with a resolution of 0.649. Based on the  $1/K_0$  values, our results indicate that the PG has a more compact structure than the BMP, which is consistent with previous studies.<sup>34,59</sup>

#### Acyl Chain Composition, Double Bond Position and Double Bond Geometry


The biophysical effects of acyl chain compositions and the positions and geometries (*cis* (Z) vs *trans* (E)) of double bonds in lipid acyl chains are well documented.<sup>60,61</sup> The thickness and fluidity of biological membranes are affected by lipid double bond geometry, as well as acyl chain length.<sup>32,34</sup> Specifically, shorter acyl chains and *cis*- double bonds result in thinner, more fluidic membranes.<sup>62</sup> Lipids with *trans* double bonds are known to cause changes in membrane function and inhibition of lipid enzymatic pathways<sup>32</sup>. Differences in double bond position, on the other hand, have been linked to the development of various

diseases, including cancer and type 2 diabetes,<sup>63</sup> and have been known to influence lipid-protein binding.<sup>64</sup> To assess the separation capabilities of MALDI TIMS for isomers in both negative and positive ionization mode, lipids with differences in the acyl chain composition, double bond position, double bond geometry, and combinations of composition, position, and geometry were examined. First, MALDI TIMS was performed in negative ionization mode for the deprotonated ions of PG (18:1(9E)/18:1(9E)) and PG (18:0/18:2(9Z,12Z)) (Figure 3B). Both lipids have two unsaturations: PG (18:0/18:2(9Z,12Z)) has two double bonds on the same acyl chain in a *cis* conformation and PG (18:1/18:1(9E)) has a single unsaturation on each chain in the *trans* geometry; these lipid standards differ in acyl chain composition, double bond position, and double bond geometry. The mass spectrum of each standard was dominated by the  $[M-H]^-$  ion at  $m/z$  773.53. The measured  $1/K_0$  values for PG (18:1(9E)/18:1(9E)) and PG (18:0/18:2(9Z,12Z)) were 1.457 and 1.450, respectively. The lower  $1/K_0$  value for PG (18:0/18:2(9Z,12Z)) suggests it assumes a more compact gas-phase structure. A mixture of the two standards shows partial separation ( $r = 0.66$ ), suggesting that isomers with differences in acyl chain composition, double bond position, and geometry could be resolved *in situ*. To further test the separation capabilities for more subtle differences in structures, mobility data

were collected for PG (18:1(9E)/18:1(9E)) and PG (18:1(9Z)/18:1(9Z)), where the only difference is double bond geometry. Extracted ion mobilograms of the deprotonated species ( $m/z$  773.53) are shown in **Figure 3C**. The measured  $1/K_0$  values for individual PG (18:1(9Z)/18:1(9Z)) and PG (18:1(9E)/18:1(9E)) are 1.449 and 1.457, respectively. The mixture of PG (18:1(9Z)/18:1(9Z)) and PG (18:1(9Z)/18:1(9Z)) confirms that the deprotonated species of each lipid can be well resolved ( $r = 0.93$ ). Based on the lower mobility value, PG (18:1(9Z)/18:1(9Z)) takes on a more compact gas-phase structure than PG (18:1(9E)/18:1(9E)). Our results from both isomer pairs suggest that the lipids with *cis* double bonds have more compact structures than lipids with *trans* double bonds. These findings are aligned with existing literature, as *cis* double bonds are believed to introduce a bend or curl in otherwise linear acyl chains.<sup>31</sup>

MALDI TIMS separations of double bond positional and geometric isomers were also evaluated for cationic lipids. All positive ionization mode analysis was carried out on sodium-adducted PC standards, as cationization with alkali metals, has been shown to enhance ion mobility separations of isomeric species.<sup>33</sup> **Figure 3D** shows the extracted ion mobilograms of the  $[M+Na]^+$  ion ( $m/z$  808.57) for PC (18:1(6Z)/18:1(6Z)) and PC (18:1(9E)/18:1(9E)) that vary in both double bond position and double bond geometry. The individual mobilograms of the sodiated species show that PC (18:1(9E)/18:1(9E)) ( $1/K_0$  1.474) has a lower mobility value than PC (18:1(6Z)/18:1(6Z)) ( $1/K_0$  1.488); the two standards are well resolved from a mixture of the standards ( $r = 0.972$ ). To investigate the effect of double bond position alone, standard PC (18:1(6Z)/18:1(6Z)) was analyzed with PC (18:1(9Z)/18:1(9Z)) (**Figure 3E**). The individual extracted ion mobilograms ( $m/z$  808.57) show that the 9Z configuration results in a smaller structure than the 6Z, with  $1/K_0$  1.481 and 1.488, respectively. The mobilograms of the mixture shows separation of the two species with a resolution of 0.664. Our results indicate that the 9Z geometry results in a more compact structure. This finding is also supported by previous literature stating that *cis* double bonds further from the carboxyl end curl the acyl chain more, causing smaller conformations.<sup>34</sup> In addition to the aforementioned lipid isomers, PC (18:1(9E)/18:1(9E)) and PC (18:1(9Z)/18:1(9Z)) were analyzed to evaluate the separation of isomers with *cis* vs *trans* geometry (**Figure S4**). Our results suggest that the 9E configuration results in a more compact structure than the 6Z, which is consistent with previously published data suggesting that the *trans* geometry is more impacted by the sodium adduct, pulling it closer to headgroup, resulting in a more compact conformation.<sup>49</sup>

Lastly, mobility separations were performed on the  $[M+Na]^+$  ion ( $m/z$  782.56) of PC (18:1(9Z)/16:0) and PC (18:1(11Z)/16:0) (**Figure 3F**), isomers that vary only in the position of a single double bond. These lipid species are generally abundant in biological tissues and represent an interesting biological target. For instance, Paine et al. have demonstrated that although both PC (18:1(9Z)/16:0) and PC (18:1(11Z)/16:0) are found in rat brain tissue, their relative abundances vary between white matter and gray matter.<sup>65</sup> The individual mobilograms suggest that the 11Z ( $1/K_0$  1.521) isomer has a more compact structure than the 9Z (1.528). Despite the difference in a single double bond location being very subtle, slight separation in the mobilogram was still observed ( $r = 0.413$ ).

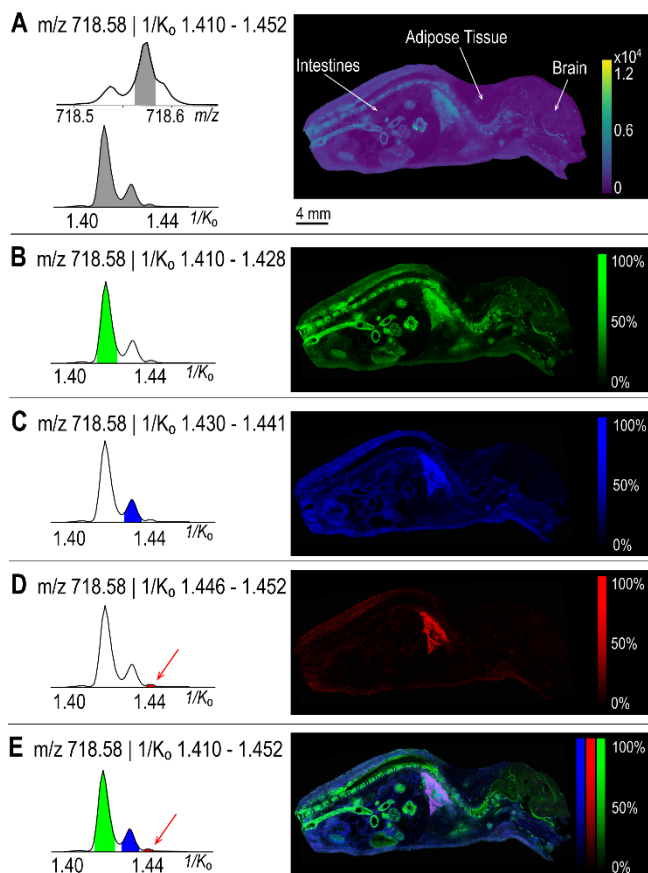

**Figure 4:** TIMS enables the separation of lipid isomers in whole-body mouse pup tissues. Composite image of all three peaks in the extracted ion mobilogram ( $1/K_0$  1.410 – 1.452) of  $m/z$  718.58 is shown (**A**). Each panel highlights an ion image generated for three distinct mobility range –  $1/K_0$  1.410-1.428 (**B**), 1.430 - 1.441 (**C**), and 1.446 – 1.452 (**D**), and an overlay of all three mobility peaks (**E**). The three ions have different localizations across the tissue sample, however, the average ion image (**A**) is dominated by the spatial distributions of the highest intensity ion (**B**). The ions were putatively identified as  $[CerP(t40:1) + H]^+$ ,  $[PC(O-32:1) + H]^+$  and  $[PC(P-32:0) + H]^+$  (0.14 ppm error) based on accurate mass alone.

### MALDI TIMS Separation and Localization of Lipid Isomers in Whole-body Mouse Pup Tissue

MALDI TIMS-IMS was used to visualize the spatial distribution of lipid isomers in whole-body mouse pup tissue. From the average mass spectrum, an example of isomeric overlap is highlighted at  $m/z$  718.58, where three distinct species are resolved in the extracted ion mobilogram (**Figure 4A**). An composite ion image for  $m/z$  718.58 including all observed mobility peaks from  $1/K_0$  1.410 – 1.452 (**Figure 4A**) was generated, followed by ion images for each ion mobility separated species (**Figures 4B-4D**). Each ion has a unique spatial distribution within the mouse pup tissue. The ion with  $1/K_0$  1.410 - 1.428 (**Figure 4B**) is the highest intensity ion that localizes to the organs of the abdominal cavity, the spinal cord, and the brain with lower intensity. The ion with  $1/K_0$  1.430 - 1.441 (**Figure 4C**) localizes primarily to the liver, to the connective tissue, and is completely absent in the brain. Lastly, the ion with mobility  $1/K_0$  1.446 - 1.452 (**Figure 4D**) localizes to the liver and has very low

intensity or is completely absent throughout the rest of the body. From the composite image, the highest intensity ion (1/K<sub>0</sub> 1.410 – 1.441) dominates the ion image, and the spatial distributions of the other two ions are lost. An overlay image of each of the ion mobility-resolved species is shown in **Figure 4E** highlighting the unique spatial distributions of each of these lipids. To determine if these ions were isobars or isomers, a serial section of the whole-body mouse pup was analyzed using ultra-high resolution FT-ICR MS (resolving power 600,000 at *m/z* 760). The whole-body average FT-ICR spectrum suggests that the ions separated in mobility space are, in fact, isomers (**Figure S5**). An ion image was generated for *m/z* 718.58, and the localizations match well to the MALDI TIMS composite ion image (**Figure 4A**). This data further confirms that the composite ion image of *m/z* 718.58 is biased toward the most intense ion in the mass spectrum (**Figure S5**). The FT-ICR MS data was also used to assist in identifying the ions (**Figure S5**). From the FT-ICR mass spectrum, there are no nominally interfering peaks and no overlapping isotopic distribution could be detected. Based on mass accuracy alone the three isomeric species were tentatively identified as [CerP(t40:1) + H]<sup>+</sup>, [PC(O-32:1) + H]<sup>+</sup> and [PC(P-32:0) + H]<sup>+</sup> (0.14 ppm error). This example highlights the utility of coupling TIMS with MALDI IMS.

## CONCLUSIONS

This work has demonstrated the use of a MALDI timsTOF IMS for enhanced peak capacity and lipid isomer analysis. An additional 669 features were detected by incorporating TIMS, resulting in a 271% increase in peak capacity. MALDI TIMS separation of lipid isomer standards, including sn-backbone isomers, acyl chain isomers, as well as double bond positional and geometric isomers was also demonstrated. Finally, *in situ* separation and imaging of lipid isomers with distinct spatial distributions was performed using whole-body mouse pup tissue. Gas-phase isobar and isomer separations confirm TIMS as a valuable tool that can bring new dimensions to IMS experiments. Lipid identification, however, remains a challenge, and TIMS data was supplemented with MALDI FT-ICR imaging to provide tentative lipid isomer identifications.

MS/MS can readily provide information on acyl chain composition for negative ion mode analysis, and head group information in positive ion mode. However, traditional tandem MS techniques cannot elucidate subtle structural differences, such as double bond position and geometry. Novel MS/MS methods, including low-energy collision-induced dissociation with Paternò-Büchi derivatization<sup>66</sup>, ozone-induced dissociation<sup>67</sup>, and ultraviolet photodissociation<sup>68</sup>, have been shown to enable a higher degree of lipid structural characterization. Supplementing the MALDI TIMS workflow with these advanced structural identification methods has the potential to further increase the depth of characterization of lipids *in situ*.

## ASSOCIATED CONTENT

### Supporting Information

Table S1. Experimental TIMS Parameters; Figure S1. Names, Structures, and Masses of lipid standards; Figure S2. Negative ionization mode mass spectra of standards; Figure S3. Positive ionization mode mass spectra of standards; Figure S4. Extracted ion mobilogram of PC (18:1(9Z)/18:1(9Z)) and PC (18:1(9E)/18:1(9E)); Figure S5. MALDI FT-ICR data.

## AUTHOR INFORMATION

### Corresponding Author

\* Jeffrey M. Spraggins, Jeff.Spraggins@vanderbilt.edu  
Jeffrey M. Spraggins: 0000-0001-9198-5498

Katerina V. Djambazova:0000-0002-2680-9014  
Lukasz G. Migas:0000-0002-1884-6405  
Dustin R. Klein: 0000-0001-7327-6479  
Emilio S. Rivera:0000-0003-3215-1092  
Elizabeth K. Neumann:0000-0002-6078-3321  
Raf Van de Plas: 0000-0002-2232-7130  
Richard M. Caprioli: 0000-0001-5859-3310

### Author Contributions

All authors have given approval to the final version of the manuscript.

## ACKNOWLEDGMENT

Support was provided by The National Science Foundation Major Research Instrument Program (CBET – 1828299 awarded to J.M.S. and R.M.C.) and The National Institutes of Health (NIH) National Institute of General Medical Sciences (P41 GM103391 awarded to R.M.C.). The 15T FT-ICR MS in the Mass Spectrometry Research Center at Vanderbilt University was acquired through the NIH Shared Instrumentation Grant Program (1S10OD012359). E.K.N. is supported by a National Institute of Environmental Health Sciences training grant (T32ES007028). The authors wish to thank Terry Dermody from the University of Pittsburgh for providing the mouse pup and Daniel Ryan for sample assistance.

## REFERENCES

- (1) Shevchenko, A.; Simons, K. Lipidomics: Coming to Grips with Lipid Diversity. *Nat. Rev. Mol. Cell. Biol.* **2010**, *11* (8), 593–598. <https://doi.org/10.1038/nrm2934>.
- (2) Liebisch, G.; Vizcaíno, J. A.; Köfeler, H.; Trötzmüller, M.; Griffiths, W. J.; Schmitz, G.; Spener, F.; Wakelam, M. J. O. Shorthand Notation for Lipid Structures Derived from Mass Spectrometry. *J. Lipid Res.* **2013**, *54* (6), 1523–1530. <https://doi.org/10.1194/jlr.M033506>.
- (3) Yang, K.; Han, X. Lipidomics: Techniques, Applications, and Outcomes Related to Biomedical Sciences. *Trends Biochem. Sci.* **2016**, *41* (11), 954–969. <https://doi.org/10.1016/j.tibs.2016.08.010>.
- (4) O'Donnell, V. B.; Ekroos, K.; Liebisch, G.; Wakelam, M. Lipidomics: Current State of the Art in a Fast Moving Field. *Syst. Biol. Med.* **2020**, *12* (1), e1466. <https://doi.org/10.1002/wsbm.1466>.
- (5) Rustam, Y. H.; Reid, G. E. Analytical Challenges and Recent Advances in Mass Spectrometry Based Lipidomics. *Anal. Chem.* American Chemical Society January 2, 2018, pp 374–397. <https://doi.org/10.1021/acs.analchem.7b04836>.
- (6) Caprioli, R. M.; Farmer, T. B.; Gile, J. Molecular Imaging of Biological Samples: Localization of Peptides and Proteins Using MALDI-TOF MS. *Anal. Chem.* **1997**, *69* (23), 4751–4760. <https://doi.org/10.1021/AC970888I>.
- (7) Angel, P. M.; Spraggins, J. M.; Baldwin, H. S.; Caprioli, R. Enhanced Sensitivity for High Spatial Resolution Lipid Analysis by Negative Ion Mode Matrix Assisted Laser Desorption Ionization Imaging Mass Spectrometry. *Anal. Chem.* **2012**, *84* (3), 1557–1564. <https://doi.org/10.1021/ac202383m>.
- (8) Berry, K. A. Z.; Hankin, J. A.; Barkley, R. M.; Spraggins, J. M.; Caprioli, R. M.; Murphy, R. C. MALDI Imaging of Lipid Biochemistry in Tissues by Mass Spectrometry. *Chem. Rev.* **2011**, *111* (10), 6491. <https://doi.org/10.1021/CR200280P>.
- (9) Cohen, L. H.; Gusev, A. I. Small Molecule Analysis by MALDI Mass Spectrometry. *Anal. Bioanal. Chem.* **2002**, *373* (7), 571–586. <https://doi.org/10.1007/s00216-002-1321-z>.
- (10) Bhandari, D. R. a.; Schott, M.; Römpf, A.; Vilcinskas, A.; Spengler, B. Metabolite Localization by Atmospheric Pressure



- High-Resolution Scanning Microprobe Matrix-Assisted Laser Desorption/Ionization Mass Spectrometry Imaging in Whole-Body Sections and Individual Organs of the Rove Beetle *Pederus Riparius*. *Anal. Bioanal. Chem.* **2015**, *407* (8), 2189–2201. <https://doi.org/10.1007/s00216-014-8327-1>.
- (11) Powers, T. W.; Neely, B. A.; Shao, Y.; Tang, H.; Troyer, D. A.; Mehta, A. S.; Haab, B. B.; Drake, R. R. MALDI Imaging Mass Spectrometry Profiling of N-Glycans in Formalin-Fixed Paraffin Embedded Clinical Tissue Blocks and Tissue Microarrays. *PLoS One* **2014**, *9* (9), e106255. <https://doi.org/10.1371/journal.pone.0106255>.
- (12) Drake, R. R.; Powers, T. W.; Norris-Caneda, K.; Mehta, A. S.; Angel, P. M. In Situ Imaging of N-Glycans by MALDI Imaging Mass Spectrometry of Fresh or Formalin-Fixed Paraffin-Embedded Tissue. *Curr. Protoc. Protein Sci.* **2018**, *94* (1), e68. <https://doi.org/10.1002/cpps.68>.
- (13) Schober, Y.; Guenther, S.; Spengler, B.; Römpf, A. High-Resolution Matrix-Assisted Laser Desorption/Ionization Imaging of Tryptic Peptides from Tissue. *Rapid Commun. Mass Spectrom.* **2012**, *26* (9), 1141–1146. <https://doi.org/10.1002/rcm.6192>.
- (14) Prentice, B. M.; Ryan, D. J.; Van De Plas, R.; Caprioli, R. M.; Spraggins, J. M. Enhanced Ion Transmission Efficiency up to  $m/z$  24 000 for MALDI Protein Imaging Mass Spectrometry. *Anal. Chem.* **2018**, *90* (8), 5090–5099. <https://doi.org/10.1021/acs.analchem.7b05105>.
- (15) Spraggins, J. M.; Caprioli, R. M. High-Speed MALDI-TOF Imaging Mass Spectrometry: Rapid Ion Image Acquisition and Considerations for Next Generation Instrumentation. *J. Am. Soc. Mass Spectrom.* **2011**, *22* (6), 1022–1031. <https://doi.org/10.1007/s13361-011-0121-0>.
- (16) El-Aneel, A.; Cohen, A.; Banoub, J. Mass Spectrometry, Review of the Basics: Electrospray, MALDI, and Commonly Used Mass Analyzers. *Appl. Spectrosc. Rev.* **2009**, *44* (3), 210–230. <https://doi.org/10.1080/05704920902717872>.
- (17) Smith, D. F.; Kiss, A.; Leach, F. E.; Robinson, E. W.; Paša-Tolić, L.; Heeren, R. M. A. High Mass Accuracy and High Mass Resolving Power FT-ICR Secondary Ion Mass Spectrometry for Biological Tissue Imaging. *Anal. Bioanal. Chem.* **2013**, *405* (18), 6069–6076. <https://doi.org/10.1007/s00216-013-7048-1>.
- (18) Kompauer, M.; Heiles, S.; Spengler, B. Atmospheric Pressure MALDI Mass Spectrometry Imaging of Tissues and Cells at 1.4-Mm Lateral Resolution. *Nat. Methods* **2017**, *14* (1), 90–96. <https://doi.org/10.1038/nmeth.4071>.
- (19) Porta Siegel, T.; Ekroos, K.; Ellis, S. R. Reshaping Lipid Biochemistry by Pushing Barriers in Structural Lipidomics. *Angew. Chem. Int. Ed.* **2019**, *58* (20), 6492–6501. <https://doi.org/10.1002/anie.201812698>.
- (20) Damen, C. W. N.; Isaac, G.; Langridge, J.; Hankemeier, T.; Vreeken, R. J. Enhanced Lipid Isomer Separation in Human Plasma Using Reversed-Phase UPLC with Ion-Mobility/High-Resolution MS Detection. *J. Lipid Res.* **2014**, *55* (8), 1772–1783. <https://doi.org/10.1194/jlr.D047795>.
- (21) May, J. C.; Goodwin, C. R.; Lareau, N. M.; Leaprot, K. L.; Morris, C. B.; Kurulugama, R. T.; Mordehai, A.; Klein, C.; Barry, W.; Darland, E.; et al. Conformational Ordering of Biomolecules in the Gas Phase: Nitrogen Collision Cross Sections Measured on a Prototype High Resolution Drift Tube Ion Mobility-Mass Spectrometer. *Anal. Chem.* **2014**, *86* (4), 2107–2116. <https://doi.org/10.1021/ac4038448>.
- (22) Valentine, S. J.; Koeniger, S. L.; Clemmer, D. E. A Split-Field Drift Tube for Separation and Efficient Fragmentation of Biomolecular Ions. *Anal. Chem.* **2003**, *75* (22), 6202–6208. <https://doi.org/10.1021/ac030111r>.
- (23) Zhong, Y.; Hyung, S.-J.; Ruotolo, B. T. Characterizing the Resolution and Accuracy of a Second-Generation Traveling-Wave Ion Mobility Separator for Biomolecular Ions. *Analyst* **2011**, *136* (17), 3534. <https://doi.org/10.1039/c0an00987c>.
- (24) Shvartsburg, A. A.; Smith, R. D. Fundamentals of Traveling Wave Ion Mobility Spectrometry. *Anal. Chem.* **2008**, *80* (24), 9689–9699. <https://doi.org/10.1021/ac8016295>.
- (25) Deng, L.; Ibrahim, Y. M.; Baker, E. S.; Aly, N. A.; Hamid, A. M.; Zhang, X.; Zheng, X.; Garimella, S. V. B.; Webb, I. K.; Prost, S. A.; et al. Ion Mobility Separations of Isomers Based upon Long Path Length Structures for Lossless Ion Manipulations Combined with Mass Spectrometry. *ChemistrySelect* **2016**, *1* (10), 2396–2399. <https://doi.org/10.1002/slct.201600460>.
- (26) Ibrahim, Y. M.; Hamid, A. M.; Deng, L.; Garimella, S. V. B.; Webb, I. K.; Baker, E. S.; Smith, R. D. New Frontiers for Mass Spectrometry Based upon Structures for Lossless Ion Manipulations. *Analyst* **2017**, *142* (7), 1010–1021. <https://doi.org/10.1039/c7an00031f>.
- (27) Purves, R. W.; Guevremont, R. Electrospray Ionization High-Field Asymmetric Waveform Ion Mobility Spectrometry-Mass Spectrometry. *Int. J. Mass Spectrom. Ion Process.* **1970**, *52* (3), 2346–2357. <https://doi.org/10.1021/ac981380y>.
- (28) Kolakowski, B. M.; Mester, Z. Review of Applications of High-Field Asymmetric Waveform Ion Mobility Spectrometry (FAIMS) and Differential Mobility Spectrometry (DMS). *Analyst* **2007**, *132* (9), 842. <https://doi.org/10.1039/b706039d>.
- (29) Guevremont, R. High-Field Asymmetric Waveform Ion Mobility Spectrometry: A New Tool for Mass Spectrometry. *J. Chromatogr. A.* **2004**, *1058* (1–2), 3–19. <https://doi.org/10.1016/j.chroma.2004.08.119>.
- (30) Fernandez-Lima, F.; Kaplan, D. A.; Suetering, J.; Park, M. A. Gas-Phase Separation Using a Trapped Ion Mobility Spectrometer. *Int. J. Ion Mobil. Spectrom.* **2011**, *14* (2–3). <https://doi.org/10.1007/s12127-011-0067-8>.
- (31) Michelmann, K.; Silveira, J. A.; Ridgeway, M. E.; Park, M. A. Fundamentals of Trapped Ion Mobility Spectrometry. *J. Am. Soc. Mass Spectrom.* **2015**, *26* (1), 14–24. <https://doi.org/10.1007/s13361-014-0999-4>.
- (32) Grossl, M.; Graf, S.; Knochenmuss, R. High Resolution Ion Mobility-Mass Spectrometry for Separation and Identification of Isomeric Lipids. *Analyst* **2015**, *140* (20), 6904–6911. <https://doi.org/10.1039/C5AN00838G>.
- (33) Zheng, X.; Smith, R. D.; Baker, E. S. Recent Advances in Lipid Separations and Structural Elucidation Using Mass Spectrometry Combined with Ion Mobility Spectrometry, Ion-Molecule Reactions and Fragmentation Approaches. *Curr. Opin. Chem. Biol.* Elsevier Ltd February 1, 2018, pp 111–118. <https://doi.org/10.1016/j.cbpa.2017.11.009>.
- (34) Kyle, J. E.; Zhang, X.; Weitz, K. K.; Monroe, M. E.; Ibrahim, Y. M.; Moore, R. J.; Cha, J.; Sun, X.; Lovelace, E. S.; Wagoner, J.; et al. Uncovering Biologically Significant Lipid Isomers with Liquid Chromatography, Ion Mobility Spectrometry and Mass Spectrometry. *Analyst* **2016**, *141* (5), 1649–1659. <https://doi.org/10.1039/C5AN02062J>.
- (35) McLean, J. A.; Ridenour, W. B.; Caprioli, R. M. Profiling and Imaging of Tissues by Imaging Ion Mobility-Mass Spectrometry. *J. Mass Spectrom.* **2007**, *42* (8), 1099–1105. <https://doi.org/10.1002/jms.1254>.
- (36) Škrásková, K.; Claude, E.; Jones, E. A.; Towers, M.; Ellis, S. R.; Heeren, R. M. A. Enhanced Capabilities for Imaging Gangliosides in Murine Brain with Matrix-Assisted Laser Desorption/Ionization and Desorption Electrospray Ionization Mass Spectrometry Coupled to Ion Mobility Separation. *Methods* **2016**, *104*, 69–78. <https://doi.org/10.1016/j.ymeth.2016.02.014>.
- (37) Wojcik, R.; Webb, I. K.; Deng, L.; Garimella, S. V. B.; Prost, S. A.; Ibrahim, Y. M.; Baker, E. S.; Smith, R. D. Lipid and Glycolipid Isomer Analyses Using Ultra-High Resolution Ion Mobility Spectrometry Separations. *Int. J. Mol. Sci.* **2017**, *18* (1), 183. <https://doi.org/10.3390/ijms18010183>.
- (38) Ma, X.; Chong, L.; Tian, R.; Shi, R.; Hu, T. Y.; Ouyang, Z.; Xia, Y. Identification and Quantitation of Lipid C=C Location Isomers: A Shotgun Lipidomics Approach Enabled by Photochemical Reaction. *Proc. Nat. Acad. Sci.* **2016**, *113* (10), 2573–2578. <https://doi.org/10.1073/pnas.1523356113>.
- (39) Bowman, A. P.; Abzalimov, R. R.; Shvartsburg, A. A. Broad Separation of Isomeric Lipids by High-Resolution Differential Ion Mobility Spectrometry with Tandem Mass Spectrometry. *J. Am. Soc. Mass Spectrom.* **2017**, *28* (8), 1552–1561. <https://doi.org/10.1007/s13361-017-1675-2>.
- (40) Nagy, G.; Veličković, D.; Chu, R. K.; Carrell, A. A.; Weston, D. J.; Ibrahim, Y. M.; Anderton, C. R.; Smith, R. D. Towards Resolving the Spatial Metabolome with Unambiguous Molecular Annotations in Complex Biological Systems by Coupling Mass Spectrometry Imaging with Structures for Lossless Ion Manipulations. *Chem. Commun.* **2019**, *55* (3), 306–309. <https://doi.org/10.1039/c8cc07482h>.

- (41) Sans, M.; Feider, C. L.; Eberlin, L. S. Advances in Mass Spectrometry Imaging Coupled to Ion Mobility Spectrometry for Enhanced Imaging of Biological Tissues. *Curr. Opin. Chem. Biol.* **2018**, *42*, 138–146. <https://doi.org/10.1016/j.cbpa.2017.12.005>.
- (42) Shvartsburg, A. A.; Isaac, G.; Leveque, N.; Smith, R. D.; Metz, T. O. Separation and Classification of Lipids Using Differential Ion Mobility Spectrometry. *J. Am. Soc. Mass Spectrom.* **2011**, *22* (7), 1146–1155. <https://doi.org/10.1007/s13361-011-0114-z>.
- (43) Lintonen, T. P. I.; Baker, P. R. S.; Suoniemi, M.; Ubhi, B. K.; Koistinen, K. M.; Duchoslav, E.; Campbell, J. L.; Ekroos, K. Differential Mobility Spectrometry-Driven Shotgun Lipidomics. *Anal. Chem.* **2014**, *86* (19), 9662–9669. <https://doi.org/10.1021/ac5021744>.
- (44) Feider, C. L.; Elizondo, N.; Eberlin, L. S. Ambient Ionization and FAIMS Mass Spectrometry for Enhanced Imaging of Multiply Charged Molecular Ions in Biological Tissues. *Anal. Chem.* **2016**, *88* (23), 11533–11541. <https://doi.org/10.1021/acs.analchem.6b02798>.
- (45) Spraggins, J. M.; Djambazova, K. V.; Rivera, E. S.; Migas, L. G.; Neumann, E. K.; Fuetterer, A.; Suetering, J.; Goedecke, N.; Ly, A.; Van De Plas, R.; et al. High-Performance Molecular Imaging with MALDI Trapped Ion-Mobility Time-of-Flight (TimsTOF) Mass Spectrometry. *Anal. Chem.* **2019**, *91* (22), 14552–14560. <https://doi.org/10.1021/acs.analchem.9b03612>.
- (46) Ridgeway, M. E.; Lubeck, M.; Jordens, J.; Mann, M.; Park, M. A. Trapped Ion Mobility Spectrometry: A Short Review. *Int. J. Mass Spectrom.* **2018**, *425*, 22–35. <https://doi.org/10.1016/j.ijms.2018.01.006>.
- (47) Hernandez, D. R.; DeBord, J. D.; Ridgeway, M. E.; Kaplan, D. A.; Park, M. A.; Fernandez-Lima, F. Ion Dynamics in a Trapped Ion Mobility Spectrometer. *Analyst.* **2014**, *139* (8), 1913–1921. <https://doi.org/10.1039/C3AN02174B>.
- (48) Silveira, J. A.; Ridgeway, M. E.; Park, M. A. High Resolution Trapped Ion Mobility Spectrometry of Peptides. *Anal. Chem.* **2014**, *86* (12), 5624–5627. <https://doi.org/10.1021/ac501261h>.
- (49) Fouque, K. J. D.; Ramirez, C. E.; Lewis, R. L.; Koelmel, J. P.; Garrett, T. J.; Yost, R. A.; Fernandez-Lima, F. Effective Liquid Chromatography-Trapped Ion Mobility Spectrometry-Mass Spectrometry Separation of Isomeric Lipid Species. *Anal. Chem.* **2019**, *91* (8), 5021–5027. <https://doi.org/10.1021/acs.analchem.8b04979>.
- (50) Fu, T.; Oetjen, J.; Chapelle, M.; Verdu, A.; Szesny, M.; Chaumot, A.; Degli-Esposti, D.; Geffard, O.; Clément, Y.; Salvador, A.; et al. In Situ Isobaric Lipid Mapping by MALDI-Ion Mobility Separation-Mass Spectrometry Imaging. *J. Mass Spectrom.* **2020**, e4531. <https://doi.org/10.1002/jms.4531>.
- (51) Fahy, E.; Sud, M.; Cotter, D.; Subramaniam, S. LIPID MAPS Online Tools for Lipid Research. *Nucleic Acids Res.* **2007**, *35*. <https://doi.org/10.1093/nar/gkm324>.
- (52) Fahy, E.; Subramaniam, S.; Murphy, R. C.; Nishijima, M.; Raetz, C. R. H.; Shimizu, T.; Spener, F.; Van Meer, G.; Wakelam, M. J. O.; Dennis, E. A. Update of the LIPID MAPS Comprehensive Classification System for Lipids. *J. Lipid Res.* **2009**, *50* (SUPPL.), S9–S14. <https://doi.org/10.1194/jlr.R800095-JLR200>.
- (53) Bielow, C.; Mastrobuoni, G.; Orioli, M.; Kempa, S. On Mass Ambiguities in High-Resolution Shotgun Lipidomics. *Anal. Chem.* **2017**, *89* (5), 2986–2994. <https://doi.org/10.1021/acs.analchem.6b04456>.
- (54) Wang, M.; Huang, Y.; Han, X. Accurate Mass Searching of Individual Lipid Species Candidates from High-Resolution Mass Spectra for Shotgun Lipidomics. *Rapid Commun. Mass Spectrom.* **2014**, *28* (20), 2201–2210. <https://doi.org/10.1002/rcm.7015>.
- (55) Jeanne Dit Fouque, K.; Ramirez, C. E.; Lewis, R. L.; Koelmel, J. P.; Garrett, T. J.; Yost, R. A.; Fernandez-Lima, F. Effective Liquid Chromatography-Trapped Ion Mobility Spectrometry-Mass Spectrometry Separation of Isomeric Lipid Species. *Anal. Chem.* **2019**, *91* (8), 5021–5027. <https://doi.org/10.1021/acs.analchem.8b04979>.
- (56) Anderson, D. M. G.; Ablonczy, Z.; Koutalos, Y.; Hanneken, A. M.; Spraggins, J. M.; Calcutt, M. W.; Crouch, R. K.; Caprioli, R. M.; Schey, K. L. Bis(Monoacylglycerol)Phosphate Lipids in the Retinal Pigment Epithelium Implicate Lysosomal/Endosomal Dysfunction in a Model of Stargardt Disease and Human Retinas. *Sci. Rep.* **2017**, *7* (1), 17352. <https://doi.org/10.1038/s41598-017-17402-1>.
- (57) Akgoc, Z.; Sena-Esteves, M.; Martin, D. R.; Han, X.; D’Azzo, A.; Seyfried, T. N. Bis(Monoacylglycerol)Phosphate: A Secondary Storage Lipid in the Gangliosidoses. *J. Lipid Res.* **2015**, *56* (5), 1006–1013. <https://doi.org/10.1194/jlr.M057851>.
- (58) Numata, M.; Chu, H. W.; Dakhama, A.; Voelker, D. R. Pulmonary Surfactant Phosphatidylglycerol Inhibits Respiratory Syncytial Virus-Induced Inflammation and Infection. *Proc. Nat. Acad. Sci.* **2010**, *107* (1), 320–325. <https://doi.org/10.1073/pnas.0909361107>.
- (59) Hankin, J. A.; Murphy, R. C.; Barkley, R. M.; Gijón, M. A. Ion Mobility and Tandem Mass Spectrometry of Phosphatidylglycerol and Bis(Monoacylglycerol) Phosphate (BMP). *Int. J. Mass Spectrom.* **2015**, *378*, 255–263. <https://doi.org/10.1016/j.ijms.2014.08.026>.
- (60) Janmey, P. A.; Kinnunen, P. K. J. Biophysical Properties of Lipids and Dynamic Membranes. *Trends Cell. Biol.* **2006**, *16* (10), 538–546. <https://doi.org/10.1016/j.tcb.2006.08.009>.
- (61) Sandra, K.; Sandra, P. Lipidomics from an Analytical Perspective. *Curr. Opin. Chem. Biol.* **2013**, *17* (5), 847–853. <https://doi.org/10.1016/j.cbpa.2013.06.010>.
- (62) Lewis, B. A.; Engelman, D. M. Lipid Bilayer Thickness Varies Linearly with Acyl Chain Length in Fluid Phosphatidylcholine Vesicles. *J. Mol. Biol.* **1983**, *166* (2), 211–217. [https://doi.org/10.1016/S0022-2836\(83\)80007-2](https://doi.org/10.1016/S0022-2836(83)80007-2).
- (63) Mitchell, T. W.; Pham, H.; Thomas, M. C.; Blanksby, S. J. Identification of Double Bond Position in Lipids: From GC to OzID. *J. Chromatogr. B.* **2009**, *877* (26), 2722–2735. <https://doi.org/10.1016/j.jchromb.2009.01.017>.
- (64) Brown, S. H. J.; Mitchell, T. W.; Oakley, A. J.; Pham, H. T.; Blanksby, S. J. Time to Face the Fats: What Can Mass Spectrometry Reveal about the Structure of Lipids and Their Interactions with Proteins? Lipid-Protein Interactions. *J. Am. Soc. Mass Spectrom.* **2012**, *23*, 1441–1449. <https://doi.org/10.1007/s13361-012-0410-2>.
- (65) Paine, M. R. L.; Poad, B. L. J.; Eijkel, G. B.; Marshall, D. L.; Blanksby, S. J.; Heeren, R. M. A.; Ellis, S. R. Mass Spectrometry Imaging with Isomeric Resolution Enabled by Ozone-Induced Dissociation. *Angew. Chem. Int. Ed.* **2018**, *57* (33), 10530–10534. <https://doi.org/10.1002/anie.201802937>.
- (66) Ma, X.; Xia, Y. Pinpointing Double Bonds in Lipids by Paternò-Büchi Reactions and Mass Spectrometry. *Angew. Chem. Int. Ed.* **2014**, *53* (10), 2592–2596. <https://doi.org/10.1002/anie.201310699>.
- (67) Pham, H. T.; Maccarone, A. T.; Thomas, M. C.; Campbell, J. L.; Mitchell, T. W.; Blanksby, S. J. Structural Characterization of Glycerophospholipids by Combinations of Ozone- and Collision-Induced Dissociation Mass Spectrometry: The next Step towards “Top-down” Lipidomics. *Analyst.* **2014**, *139* (1), 204–214. <https://doi.org/10.1039/C3AN01712E>.
- (68) Williams, P. E.; Klein, D. R.; Greer, S. M.; Brodbelt, J. S. Pinpointing Double Bond and Sn-Positions in Glycerophospholipids via Hybrid 193 Nm Ultraviolet Photodissociation (UVPD) Mass Spectrometry. *J. Am. Chem. Soc.* **2017**, *139* (44), 15681–15690. <https://doi.org/10.1021/jacs.7b06416>.



0612\_Manuscript\_v2.pdf (0.96 MiB)

[view on ChemRxiv](#) • [download file](#)

---

# Resolving the Complexity of Spatial Lipidomics Using MALDI TIMS Imaging Mass Spectrometry

Katerina V. Djambazova<sup>1,2</sup>, Dustin R. Klein<sup>2,3</sup>, Lukasz G. Migas<sup>4</sup>, Elizabeth K. Neumann<sup>2,3</sup>, Emilio S. Rivera<sup>2,3</sup>, Raf Van de Plas<sup>4</sup>, Richard M. Caprioli<sup>1-3,5,6</sup> and Jeffrey M. Spraggins<sup>1-3\*</sup>.

<sup>1</sup>Department of Chemistry, Vanderbilt University, 7330 Stevenson Center, Station B 351822, Nashville, TN 37235, USA

<sup>2</sup>Mass Spectrometry Research Center, Vanderbilt University, 465 21<sup>st</sup> Ave S #9160, Nashville, TN 37235, USA

<sup>3</sup>Department of Biochemistry, Vanderbilt University, 607 Light Hall, Nashville, TN 37205, USA

<sup>4</sup>Delft Center for Systems and Control, Delft University of Technology, 2628 CD Delft, The Netherlands

<sup>5</sup>Department of Pharmacology, Vanderbilt University, 2220 Pierce Avenue, Nashville, TN 37232, USA

<sup>6</sup>Department of Medicine, Vanderbilt University, 465 21<sup>st</sup> Ave S #9160, Nashville, TN 37235, USA

CORRESPONDING AUTHOR: Jeffrey M. Spraggins: [jeff.spraggins@vanderbilt.edu](mailto:jeff.spraggins@vanderbilt.edu)

---

## SUPPORTING INFORMATION

### Table of Contents

Table S1	TIMS Parameter Optimization.....	Page S-2
Figure S1	Lipid Standard Structures.....	Page S-3
Figure S2	Mass Spectra of Lipid Standards (Negative ionization mode).....	Page S-4
Figure S3	Mass Spectra of Lipid Standards (Positive ionization mode).....	Page S-5
Figure S4	Extracted Ion Mobilogram PC(18:1(9Z)/18:1(9Z)) and PC(18:1(9E)/18:1(9E)).....	Page S-6
Figure S5	MALDI FT-ICR IMS Data.....	Page S-7


**Table S1.** TIMS parameters optimized for lipid standard analysis.

<b>Lipid</b>	<b>Ion</b>	<b>t<sub>trap</sub> (ms)</b>	<b>1/K<sub>0</sub></b>	<b>V<sub>ramp</sub> (V)</b>	<b>ΔV (V)</b>	<b>Sr (V.ms<sup>-1</sup>)</b>
<b>Negative ion mode</b>						
PG (14:0/14:0)	[M-H] <sup>-</sup>	500	1.28 – 1.32	177.3 – 170.5	7	0.01
BMP (14:0/14:0)						
PG (18:0/18:2 (9Z,12Z))	[M-H] <sup>-</sup>	650	1.43 – 1.48	204.5 – 196.0	8.5	0.01
PG (18:1/18:1 (9E))						
PG (18:1/18:1 (9Z))	[M-H] <sup>-</sup>	650	1.43 – 1.48	204.5 – 196.0	8.5	0.01
PG (18:1/18:1 (9E))						
<b>Positive Ion Mode</b>						
PC (18:1/18:1 (6Z))	[M+Na] <sup>+</sup>	650	1.44 – 1.52	176.9 – 165.1	11.8	0.02
PC (18:1/18:1 (9E))						
PC (18:1/18:1 (6Z))	[M+Na] <sup>+</sup>	650	1.44 – 1.52	176.9 – 165.1	11.8	0.02
PC (18:1/18:1 (9Z))						
PC (18:1(9Z)/16:0)	[M+Na] <sup>+</sup>	650	1.52 – 1.56	182.8 – 176.9	5.9	0.01
PC (18:1(11Z)/16:0)						


## sn-backbone isomers

### PG/BMP (28:0) [M-H]<sup>-</sup> m/z 665.49

BMP 14:0/14:0

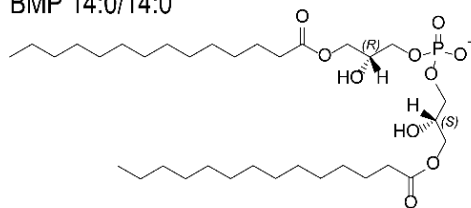


PG 14:0/14:0

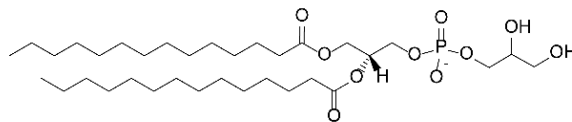


## Double bond position and geometry isomers

### PG (36:2), [M-H]<sup>-</sup> m/z 773.53


PG 18:1(9E)/18:1(9E)

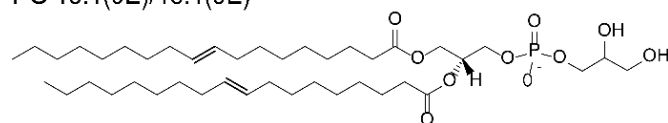

PG 18:0/18:2(9Z, 12Z)

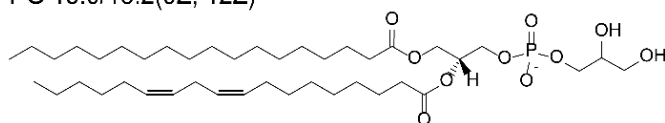


PG 18:1(9Z)/18:1(9Z)

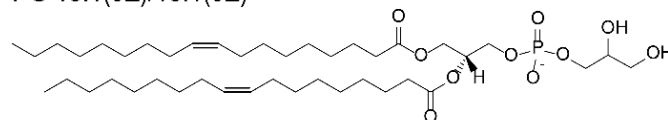


### PG (36:2), [M+Na]<sup>+</sup> m/z 808.57


PC 18:1(6Z)/18:1(6Z)

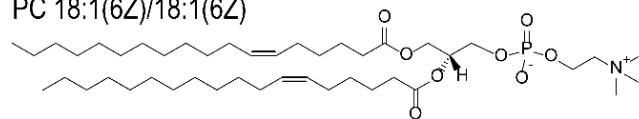

PC 18:1(9E)/18:1(9E)

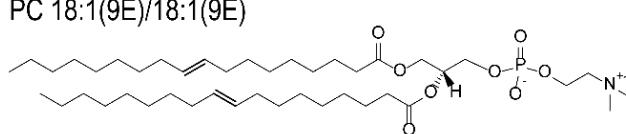


PC 18:1(9Z)/18:1(9Z)

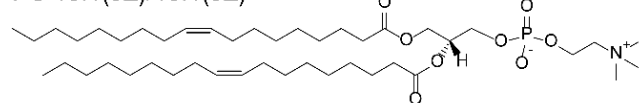


### PC (34:1), [M+Na]<sup>+</sup> m/z 782.56


PC 18:1(9Z)/16:0

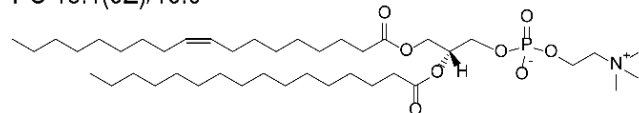

PC 18:1(11Z)/16:0

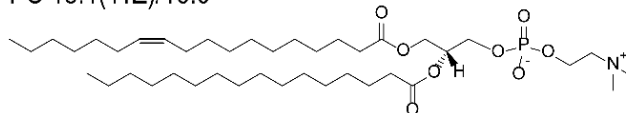


**Figure S1.** Molecular structures of lipid isomer standards analyzed by MALDI TIMS-MS

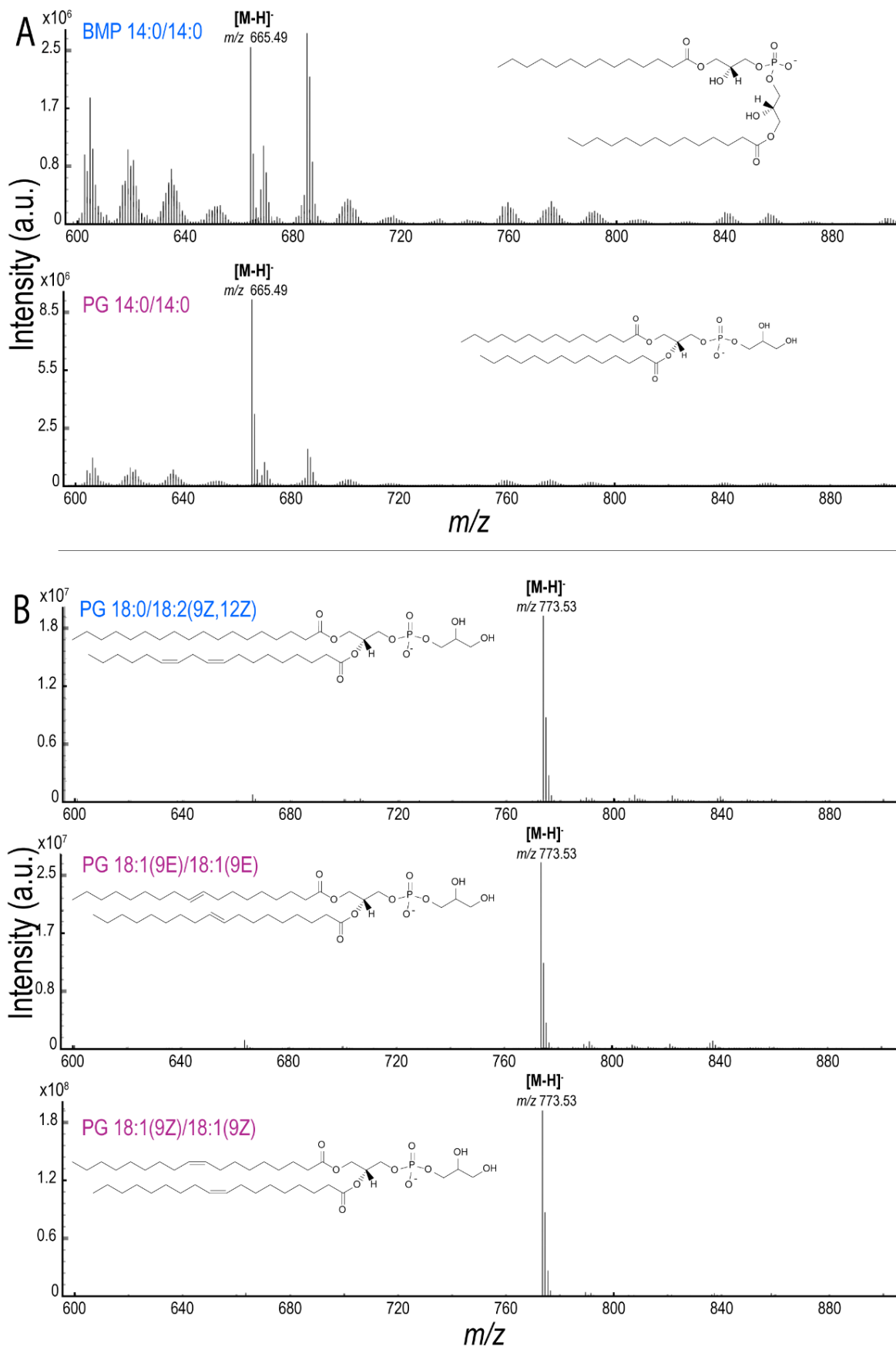


**Figure S2.** Negative ionization mode MALDI TIMS mass spectra of BMP(14:0/14:0) and PG(14:0/14:0) (**A**); PG (18:0/18:2(9Z,12Z)), PG (18:1(9E)/18:1(9E)), and PG(18:1(9Z)/18:1(9Z)) (**B**).




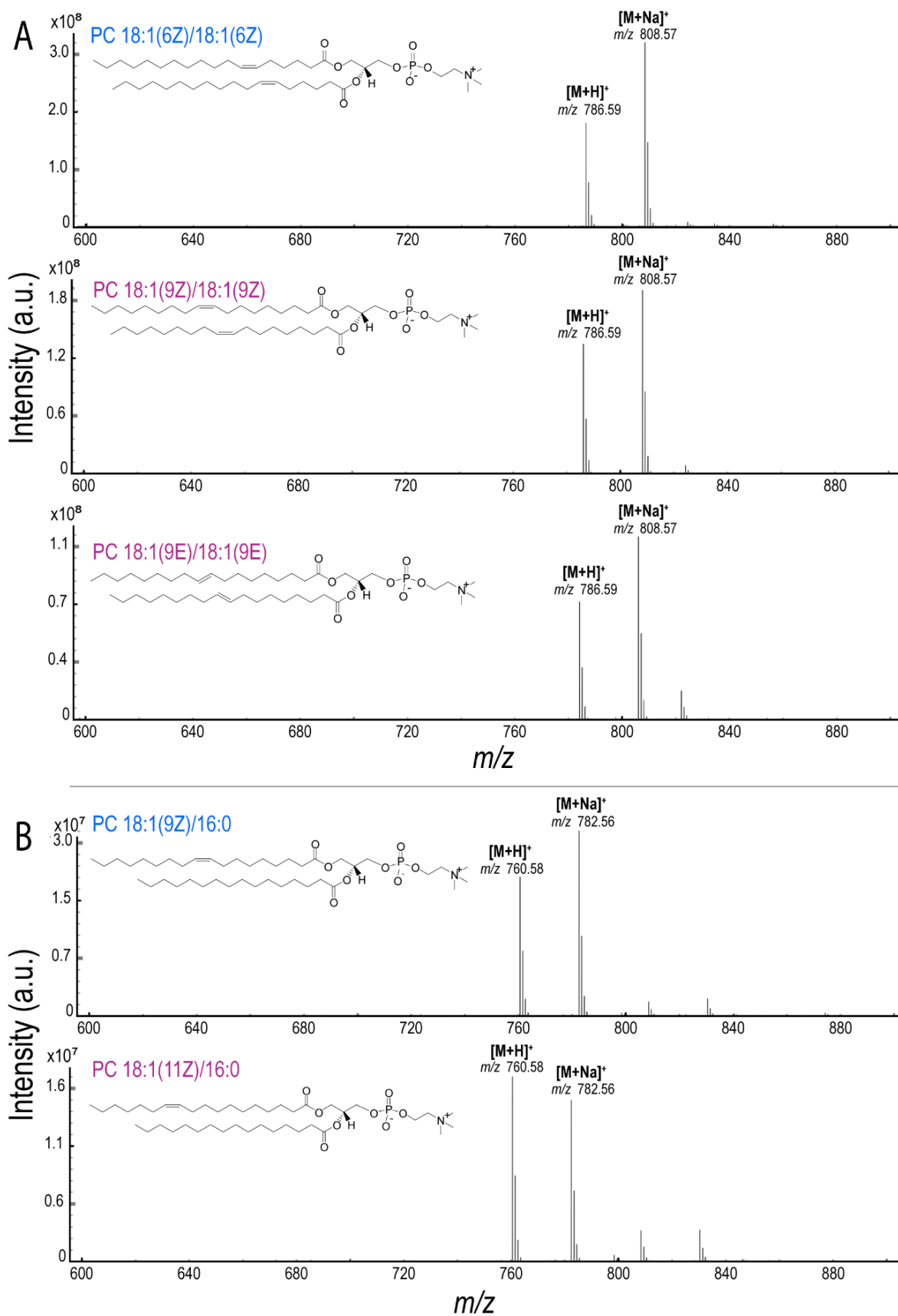

**Figure S3.** Positive ionization mode MALDI TIMS mass spectra of PC(18:1(6Z)/18:1(6Z)), PC(18:1(9Z)/18:1(9Z)), and PC(18:1(9E)/18:1(9E)) (A); PC(18:1(9Z)/16:0) and PC(18:1(11Z)/16:0) (B).

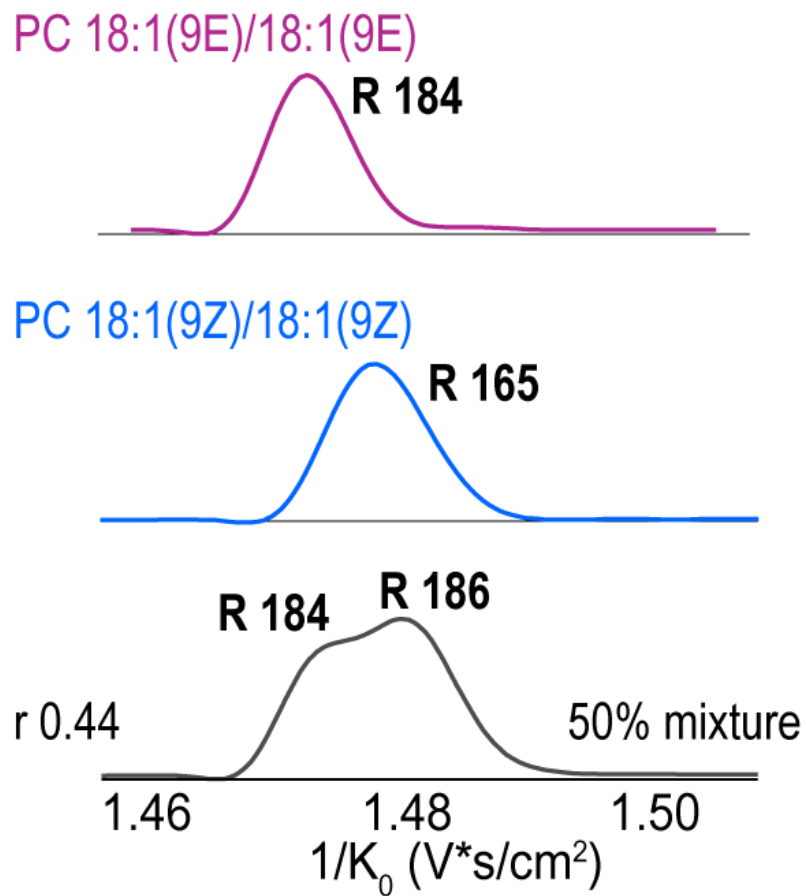

**Figure S4.** Extracted ion mobilograms of MALDI-generated ions of isomeric lipid standards for PC(18:1(9Z)/18:1(9Z)) and PC(18:1(9E)/18:1(9E)) ( $[M+Na]^+$ ,  $m/z$  808.57).

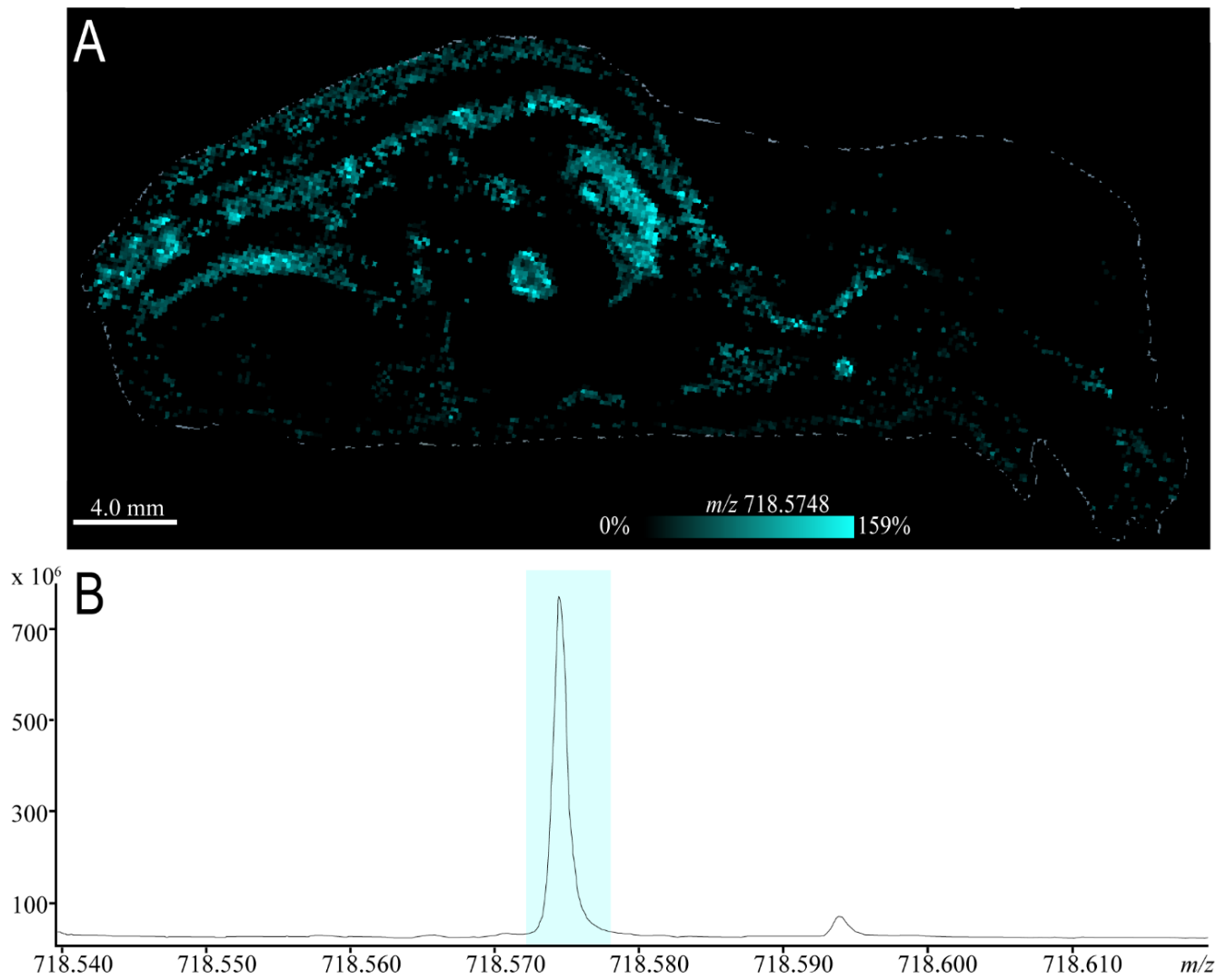

**Figure S5.** MALDI FT-ICR data generated for whole-body mouse pup. Ion image of  $m/z$  718.5748 (**A**) shows spatial localizations that mimic those available from the composite image in Figure 4A. The mass spectrum (**B**) highlights a single peak at  $m/z$  718.5748, suggesting there are no interfering peaks, due to isobaric species. The lipid isomers with  $m/z$  718.5748 were identified as  $[\text{CerP}(t40:1) + \text{H}]^+$ ,  $[\text{PC}(\text{O}-32:1) + \text{H}]^+$  and  $[\text{PC}(\text{P}-32:0) + \text{H}]^+$ .

0612\_Supplemental.pdf (889.59 KiB)

[view on ChemRxiv](#) • [download file](#)

---